




### Reproducibility (Chaney) Adapter

For consistent, reproducible injections. It also prevents the plunger from bending. The adjustable stop prevents plunger damage. Easily installed. Available for the 700, 800, 900, 1700, 1800 and 7000 series.

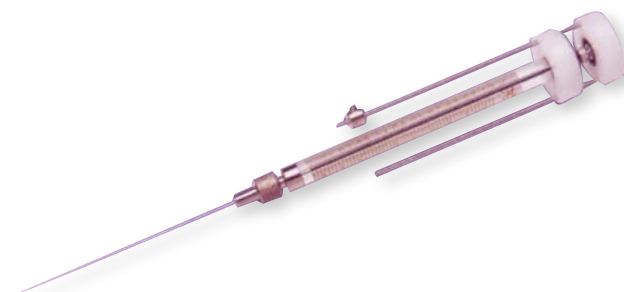

### Syringe Guide

Prevents the plunger from bending. It also prevents the plunger from being pulled out. Easily installed. Available for the 700, 1700 and 7000 series.



### MTB Sampling Valve

This lever-actuated on/off valve is used for storing and transporting samples in any GASTIGHT TLL syringe.



### Septum Adapter

Store and sample syringe contents with this adapter. It attaches to any TLL syringe or luer lock fitting. The septum is replaceable.



### Removable (RN) Adapters

There are three different removable (RN) adapters for 250  $\mu$ L and larger removable needle (RN) syringes. The port I.D. of .039" allows sample spiking with a 701N syringe. It provides connection to threaded fittings, tubing assemblies, tees, and crosses.